
# Fairywill

## Smart Oral Irrigator (F5020E) Instruction Manual

Thanks for choosing Fairywill Oral irrigator. Before operating the product, please read the manual carefully and keep it for future reference.
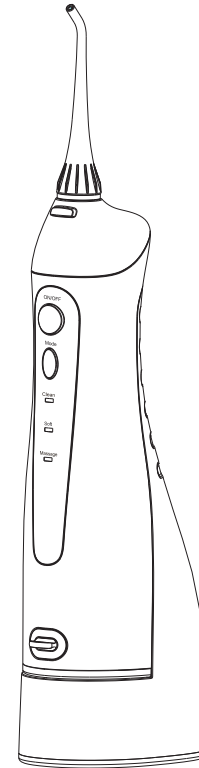
# IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

## READ ALL INSTRUCTIONS BEFORE USING (THIS APPLIANCE)

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. This appliance contains batteries that are non-replaceable. When the battery is at the end of life, the appliance shall be properly disposed of.
4. Do not expose the appliance or battery to excessive temperatures.
5. Be aware of the risk of terminals of the battery-operated appliance or battery being short circuited by metal objects.
6. The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
7. The power supply with output power not exceeding 15W in normal and abnormal working conditions must be used.
8. Toddlers, infants, people who are unable to operate the device should not use it by themselves. Please also keep it away and out of reach from children.
9. Never try to alter, disassemble or repair the device.
10. Do not fill the tank with water warmer than 40 °C.
11. If the external flexible cable or cord of this appliance is damaged, it shall be replaced by the manufacturer, a service agent or a similarly qualified person in order to avoid any hazard.
12. Keep the cord away from heated surfaces.
13. Please do not use it if you lose the nozzle.
14. Do not use for any purpose other than oral cleaning.
15. Keep nozzle away from young children. The nozzle is small enough to swallow and can be lodged in the throat.
16. Certain types of mouthwash may damage the equipment by cracking the housing or water reservoir. Therefore, do not use mouthwash in the water reservoir.
17. People who experience tooth or gum pain may not be able to use this device and should consult with their dentist before use.
18. Do not press the nozzle too firmly against your teeth or gums, otherwise it will hurt your gums.

19. Please use an adapter so that the voltage is 5V=600mA during charging, otherwise the battery could be damaged by overcharging.
20. If the battery is damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.
21. This appliance is just for household use only.

**DANGER-** To reduce the risk of electrocution while charging the appliances:

1. Do not place or store product where it can fall or be pulled into a tub or sink.
2. Do not place in or drop into water or other liquid.
3. Do not reach for a product that has fallen into water. Unplug immediately.
4. Check the charger cord for damage before the first use and during the life of the product.
5. After cleaning, make sure the charging port is completely dry before you connect it to the mains.

**WARNING -** To reduce the risk of burns, electrocution, fire, or injury to persons.

1. A product should never be left unattended when plugged in.
2. Close supervision is necessary when this product is used by, on, or near children or invalids.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.
5. Keep the cord away from heated surfaces.
6. Never use while sleeping or drowsy.
7. Never drop or insert any object into any opening or hose.
8. Do not use, store or charge the appliance next to a fire or high-temperature environment.
9. Unplug this product before filling. The reservoir with water only specified by manufacturer. Otherwise, do not overfill.


## SAVE THESE INSTRUCTIONS

There is a built-in rechargeable lithium battery in the product. When disposing of the appliance, do not mix it with general household waste. It should be sent to the battery recycling center for disposal in accordance with local disposal management regulations.

### Recycling & Environmentally Friendly Disposal

In order to protect the environment, after the expiration of the products service life, it should not be discarded in the same way as domestic garbage; please recycle at a professional recycling center. For further information on the recycling of this product, please consult your local municipality, processing service center, or the store where you purchased the product. You can help protect the environment! Please remember to respect the local regulations, and hand in the non-working electrical equipment to an appropriate waste disposal center.

## 1. Designs:

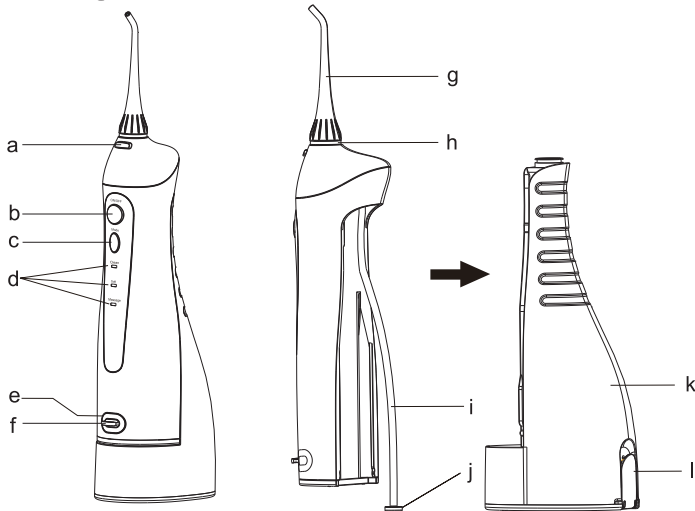- a. Nozzle release button
- b. On/Off button
- c. Mode button
- d. Mode indicator
- e. Charging base
- f. Charging base cover

- g. Nozzle
- h. Identification ring
- i. Nozzle suction tube
- j. Filter
- k. Water tank (300ML)
- l. Water tank cap

### Note:

There may be some liquid remaining inside the main device before use (tank, water suction tube, etc), this is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.

## 2. How To Use:

### 2.1 Charging the device (As Fig.1)

- Please cover the lid of the tank tightly before using the water flosser to avoid water leakage.
- Take the charger base cover down and insert the device socket (make sure no liquid is remaining inside the device plug).
- Insert the charging adaptor securely into the power outlet.
- The “CH” indicator blinks when the device is charging, and it will turn off after the device is fully charged.
- Put the cover on the charger base again after charging is finished.

### Note:

If the “CH” indicator is blinking during or after an irrigation, that means the battery needs to be recharged. A full charge will take about 4 hours.

### 2.2 Use the Oral Irrigator

- Attach the nozzle onto the unit. (Fig.2)
- Open the cap of the water tank to fill the tank with water which should not exceed 40°C, then close the cap. (Fig.3)
- Or, slide the water tank down off the unit, fill the tank with water which should not exceed 40°C, then attach the on the unit. (Fig.4)
- Press the “Mode” button to select the mode you desired. During irrigating, you can change the working mode. (Fig.5)

**Clean:** for normal cleaning, it will last for around 90s (+/-15s).

**Soft:** for gentle cleaning of sensitive gums for around 110s(+/-15s).

**Massage:** for cleaning and massaging the gums for around 100s(+/-15s).

- Hold the unit in an upright position with the nozzle tip pointing toward your teeth.
- Keep your mouth opened slightly so that the water can flow out.

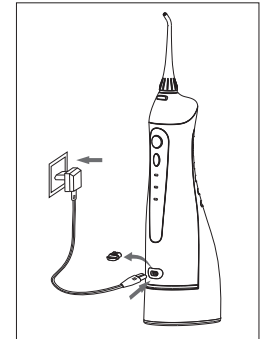
Fig.1

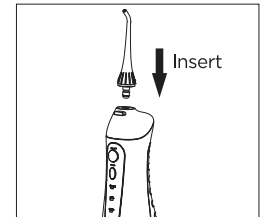
Fig.2

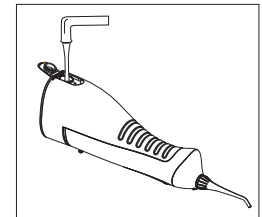
Fig.3

- Press “On/Off” button to start irrigating. (Fig.6)
- 1) Direct the water jet flow at a 90-degree angle to the teeth and gums.
- 2) Direct the water jet slowly along the teeth.
- 3) Direct the water jet along the gum line with the nozzle tip at a 90 degree angle to the gum tissue.
- 4) Cleaning between the teeth and around any braces, crowns or bridgework thoroughly.
- Press the “On/Off” button and then press the nozzle release button to replace other nozzles when needed. (Fig.7)
- Once irrigating is complete, press the “On/Off” button again to turn off the device.

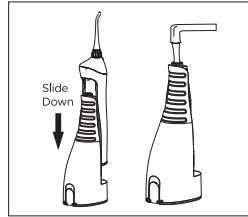
Fig.4

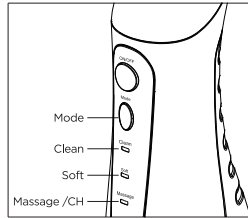
Fig.5

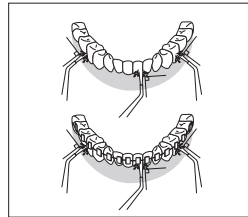
Fig.6

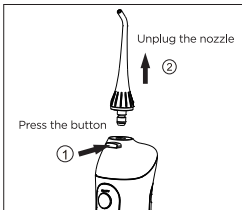
Fig.7

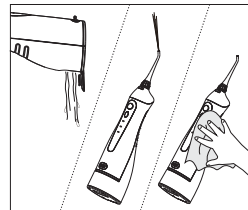
Fig.8

### 2.3 After using the appliance (Fig.8)

- Open the cap of the water tank and dispose of the water in the tank.
- Press the “On/Off” button to turn on the unit, and drain the water inside the device.
- Press the “On/Off” button again to turn off the unit, and close the water tank cap.
- Wipe the device with a dry cloth.

### Note:

To prevent bacteria from propagating, make sure there is no remaining water inside the unit. If you're not going to use this appliance for a long time, be sure to wipe and dry the appliance at your soonest convenience before storing it.

## 3. Maintenance:

The device can be only cleaned by water or a neutral cleaning agent. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may damage the device. Do not use hot water over 50°C to wash it. Keep the device away from high-temperature areas and avoid direct sunlight.

### 3.1 Main Unit

- Remove the dirt by using a cloth with liquid soap to clean the unit and a dry cloth to wipe.
- When cleaning the main unit, do not submerge it in water.
- Dry the unit after washing with water.

### 3.2 Water Tank Cap

Wipe off any grime and dirt with a cloth.

### 3.3 Nozzle

- Wash with water and wipe with a soft cloth.
- Replace the nozzle every 6 months.

### 3.4 Water Suction Tube

- Wash with water and wipe with a soft cloth.
- Do not bend, pull or twist the tube.

### 3.5 Tank

- Wash with water and wipe with soft cloth.
- Dry the tank when it is not in use for an extended period of time.

## 4. Multifunctional Nozzles Using Instruction:

**4.1 Standard Nozzles:** Suitable for daily and normal cleaning.

**4.2 Periodontal Nozzles:** Suitable for users who have sensitive teeth or gingivitis.

**4.3 Orthodontic Nozzles:** Suitable for teeth with braces or during the orthodontics period.

**4.4 Tongue Cleaning Nozzles:** Suitable for cleaning tongue coating and reducing the growth of bacteria on the surface.


## 5. Troubleshooting:

Problem	Possible Cause	Solution
It does not work	You've just purchased the product or it has been left unused for more than 3 months	Charge the unit for at least 4 hours.
It only works for a few minutes even after charging	The battery has reached the end of its lift cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the unit for at least 4 hours.
Water pressure is very low	There is dirt stuck in the nozzle	Replace the nozzle
	The nozzle is deformed	
	There is dirt stuck in the filter	Clean the filter
Water does not come out	The tank is empty	Fill up tank
	The appliance is tilted too much while in use	Use with the appliance in a vertical state

## 6. Disposal:

The rechargeable battery inside your appliance is easy to remove for recycling. At the end of the product's life and before disposal, please remove the rechargeable battery from the handle by following the instructions below. During the operation below, be sure to protect your eyes, hands, fingers, and the surface on which you work.

- Remove 2 pieces of rubber feet on the bottom of the handle (Like ① on (Fig.9))
- Insert a screwdriver into the slot located at the bottom of the handle and turn counter-clockwise to release the bottom cap. (Like ② on Fig.9)
- Slide the water tank down off the device and unscrew the two screws next to the suction tube. (Fig.10)

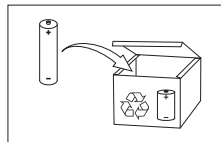- Open the top cover and put the inner part outside of the unit, then unscrew all the screws of the inner part. (Fig.11)
- Cut off the battery line and pull out the batteries from the inner part. (Fig.12)

If you have recycling questions, please contact your local waste management office.

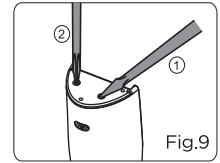
Fig.9

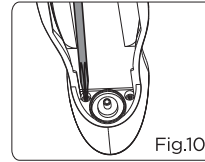
Fig.10

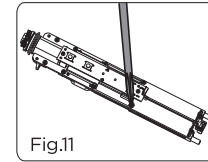
Fig.11

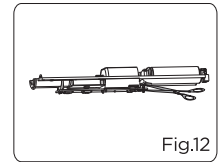
Fig.12

## 7. Warranty:

- 1-year warranty is offered, starting from the first day of purchase. You may contact Fairywill.
- During the warranty period, for any defect that may arise from materials or product process, we will replace the parts or the whole product according to certain circumstances.
- When you apply for free repair, you must submit the product and proof that can prove the date of purchase to Fairywill.
- The following conditions are not included in the warranty:
  - Damages caused by improper use, for example, wrong use of power supply, plugged into an improper socket, artificial damage etc.
  - Ageing and normal wear and tear from use.
  - Problems occur due to its use for other purposes, or improper use or use components not supplied by Fairywill.
- This warranty does not affect your statutory rights.

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