beurer

FWM 50

(D) Shiatsu-Fußwärmer Gebrauchsanweisung	2	(I) Scaldapiedi shiatsu Istruzioni per l'uso	23
(GB) Shiatsu foot warmer Instructions for use	8	(TR) Shiatsu ayak ısıtıcısı Kullanma Talimatı	28
(F) Chauffe-pieds Shiatsu Mode d'emploi	13	(RUS) Грелка для ног с массажем шиatsu Инструкция по применению.....	33
(E) Calientapiés con masaje Shiatsu Instrucciones de uso.....	18	(PL) But rozgrzewający Shiatsu Instrukcja obsługi.....	39

Inhalt

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Contents

1. About the Shiatsu massage	9	6. Cleaning and maintenance.....	12
2. Included in delivery	9	7. Storage	12
2.1 Unit description	9	8. Disposal	12
3. Important instructions	9	9. What if there are problems?	12
4. Intended use	11	10. Technical data	12
5. Operation	11		
5.1 Initial use.....	11		
5.2 Switching on/off the Shiatsu function (massage)	11		
5.3 Switching on/off the heat function (☞☞☞)....	11		
5.4 Automatic switch-off	11		

Signs and symbols

	Read instructions!		Do not wash
			Do not bleach
	Do not insert needles!		Do not tumble dry
			Do not iron
	Only use in dry rooms		Do not dry clean
IPX1	Drip proof		Must not be used by very young children (0-3 years old)
	WARNING; Warning instruction indicating a risk of injury or damage to health		
	IMPORTANT; Safety note indicating possible damage to the unit/accessory		

1. About the Shiatsu massage

Shiatsu is a type of body massage that developed in Japan over 100 years ago, based on ideas taken from traditional Chinese medicine. The basic principle is that there is a system of energy channels (meridian system) within the human body.

The objective of a Shiatsu massage is to boost the physical, emotional and mental well-being of the treated person. To this end, energy blockages in the energy channels must be removed and the self-regulating powers of the body must be stimulated.

Shiatsu therapy achieves this through the application of pressure in flowing movements along the energy channels (meridians). Although Shiatsu literally translates as 'finger pressure', treatment is not just given with the fingers but also with the heels of the hand, elbows and knees. Therapists use their whole body weight to create the necessary pressure.


Massaging a specific point on the sole of the foot can influence the corresponding organs of the body. However, we recommend studying relevant literature and asking your doctor before using point reflex zone massage to achieve the desired effect for you.

2. Included in delivery

- 1 Shiatsu foot warmer with lining
- 1 Mains part
- 1 Instructions for use

2.1 Unit description

- 1. Mains part
- 2. Cable
- 3. Connector socket
- 4. Massage button
- 5. Heat button
- 6. Shiatsu foot warmer
- 7. Lining

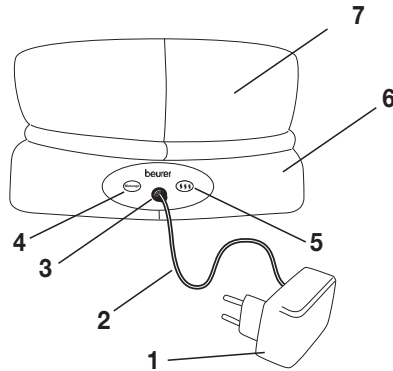3. Important instructions

Retain for future reference

WARNING

- Non-observance of the following notes may result in personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information is not only intended to protect your health and the health of others, it should also protect the product. For this reason, pay attention to these safety notes and include these instructions when handing over the product to others.
- This foot warmer must not be used by people who are not sensitive to heat or by vulnerable people who may not be able to react to overheating.
- This foot warmer must not be used by very young children (0–3 years old) as they are unable to respond to overheating.

- This foot warmer must not be used by young children (3–8 years old) unless the controller has been set by a parent or responsible adult and the child has been given adequate instructions on how to use the foot warmer safely.
- This foot warmer may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised and have been instructed on how to use the foot warmer safely, and are fully aware of the consequent risks of use.
- Children must not play with the foot warmer.
- Cleaning and user maintenance must not be performed by children unless supervised.
- This foot warmer is not designed for use in hospitals.
- This foot warmer must not be used to keep animals warm.
- Do not use this foot warmer with shoes.
- Do not insert needles.
- Only use in dry rooms. 
- Use only the original mains part provided
Type: S024FV 1200 200.
- The electrical and magnetic fields emitted by this foot warmer may interfere with the function of a pacemaker. However, they are still well within the limits: electrical field strength: max. 5000 V/m, magnetic field strength: max. 80 A/m, magnetic flux density: max. 0.1 millitesla. Please consult your doctor and the manufacturer of your pacemaker before using this foot warmer.
- Do not pull, twist or make sharp bends in the cables.
- Please check this foot warmer frequently for signs of wear and tear or damage.
If any such signs are evident, if the foot warmer has been used incorrectly or if it no longer heats up, it must be checked by the manufacturer before being switched on again.
- If the cable of the foot warmer is damaged, it must be replaced by the manufacturer, its customer service department or a similarly qualified person to avoid any risk.
- Always observe the notes on operation (section 5), cleaning and maintenance (section 6) and storage (section 7).
- If you should have any questions about using our devices, please contact Customer Services.

4. Intended use

This product is intended for warming and massaging human feet.

WARNING

Do not use the foot warmer with shoes, as they reduce your perception of heat and dirty the foot warmer. This foot warmer is not designed for use in hospitals or commercial use. The following groups in particular must not use the device: infants, small children, vulnerable persons, persons with low sensitivity to heat or animals. Otherwise, there is a risk of injury.

5. Operation

5.1 Initial use

- Connect the mains part cable to the connector socket on the back of the foot warmer.
- Insert the mains part into the wall socket.

5.2 Switching on/off the Shiatsu function (massage)

- To switch on the Shiatsu function, press the massage button on the back of the foot warmer. The massage heads will begin to rotate.
- To switch off the Shiatsu function, press the massage button on the back of the foot warmer again. The massage heads will stop.

IMPORTANT

If too much pressure is exerted on the massage heads, the massage function will remain off for safety reasons. In order to use the foot warmer again, briefly remove the cable from the connector socket and then reinsert it.

5.3 Switching on/off the heat function (☺☺☺)

- To switch on the heat function, press the heat button (☺☺☺) on the back of the foot warmer. A red LED above the heat button will light up.
- To switch off the heat function, press the heat button on the back of the foot warmer again. The red LED above the heat button will go out.

5.4 Automatic switch-off

The foot warmer will switch off after approx. 30 minutes for safety reasons. Press the massage or heat button again if you wish to continue using the foot warmer.

Switch the foot warmer off and remove the mains part from the socket if you do not wish to continue using the foot warmer after the timed switch-off.

WARNING

If the foot warmer is being used over several hours, we recommend that you select the lowest heat setting on the controller to avoid overheating the feet, which may lead to burns to the skin.

6. Cleaning and maintenance

IMPORTANT

Before cleaning, always remove the mains part from the socket first.
Otherwise, there is a risk of damage.

Small marks on the foot warmer can be removed with a damp cloth and, if necessary, with a little liquid detergent for delicate laundry.

IMPORTANT

- Please note that the foot warmer must not be washed, chemically cleaned, wrung out, tumble dried, put through a mangle or ironed. Failure to observe this may result in damage to the foot warmer.

The removable lining is machine washable at 30°.

Use a delicate laundry detergent and measure it out according to the manufacturer's instructions.

IMPORTANT

- Only use the foot warmer again once the foot warmer and lining have fully dried. Otherwise, there is a risk of damage.

7. Storage

If you do not plan to use the foot warmer for a long period of time, we recommend that you store it in the original packaging.

IMPORTANT

Please allow the foot warmer to cool down before storing. Failure to observe this may result in damage to the foot warmer. To avoid damaging the foot warmer, do not place any objects on top of it whilst storing.

8. Disposal

Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

9. What if there are problems?

Problem	Cause	Solution
The foot warmer does not warm up or the massage function cannot be activated while: – the mains part is connected to a working socket. – the mains part cable is connected to the connector socket on the foot warmer. – both functions (heat + massage) are switched on.	The foot warmer is faulty.	Send the foot warmer for servicing.

10. Technical data

See the rating label on the foot warmer.

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